

Automated Hemostasis Solution

CN-9000™


Flexible and
expandable to
optimize workflow



Together for a better
healthcare journey

Automated solutions to meet your needs

Combining high throughput and elevated efficiency, the Sysmex CN-9000 Automated Hemostasis Solution streamlines routine and specialty hemostasis testing in busy, high-volume laboratories. With continuous sample loading and automated rerun and reflex testing, the CN-9000 has a throughput of up to 1,000 samples per hour and functions as a standalone work cell or can be connected to a total laboratory automation (TLA) line.

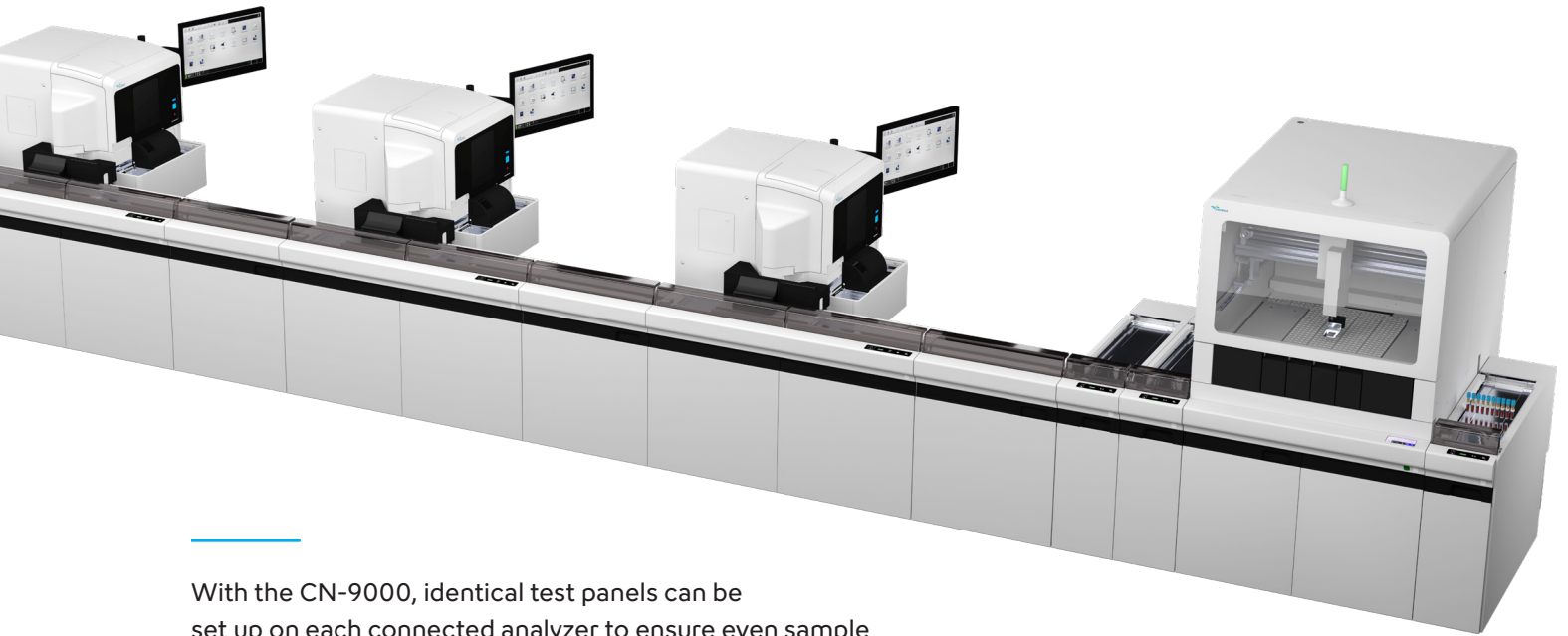

With the CN-9000, identical test panels can be set up on each connected analyzer to ensure even sample distribution on the line, or dedicated test panels can be configured on each of the analyzers to conserve reagents, QC and calibrators and labor.




Powerful productivity & exceptional service

At Sysmex, our focus is on providing solutions that make your work easier and more efficient. The CN-Series™ analyzers deliver the speed, advanced technology and reliability you need to keep pace with growing demands.

- Sysmex's fastest hemostasis solution with a throughput of up to 1000 samples per hour
- Automated repeat, re-dilute and reflex testing for walkaway operation
- Preanalytical Sample Integrity (PSI™) checks ensure sample quality
- Simultaneous multiwavelength technology for accurate results
- Reagent monitoring and robust on-board stability provide smart workflow management
- Supported by Sysmex's award-winning service teams

Save space and time

Laboratory floor space is a commodity that is not always affordable. The modular CN-9000 supports up to three CN-3000™ or CN-6000™ Blood Coagulation Analyzers as well as the optional Sysmex TS-10™ Integrated Tube Sorter and Archiver. The line can be configured with up to two turn unit modules, allowing 90° turns to the left or right, making it the 'perfect fit' for your coagulation department.

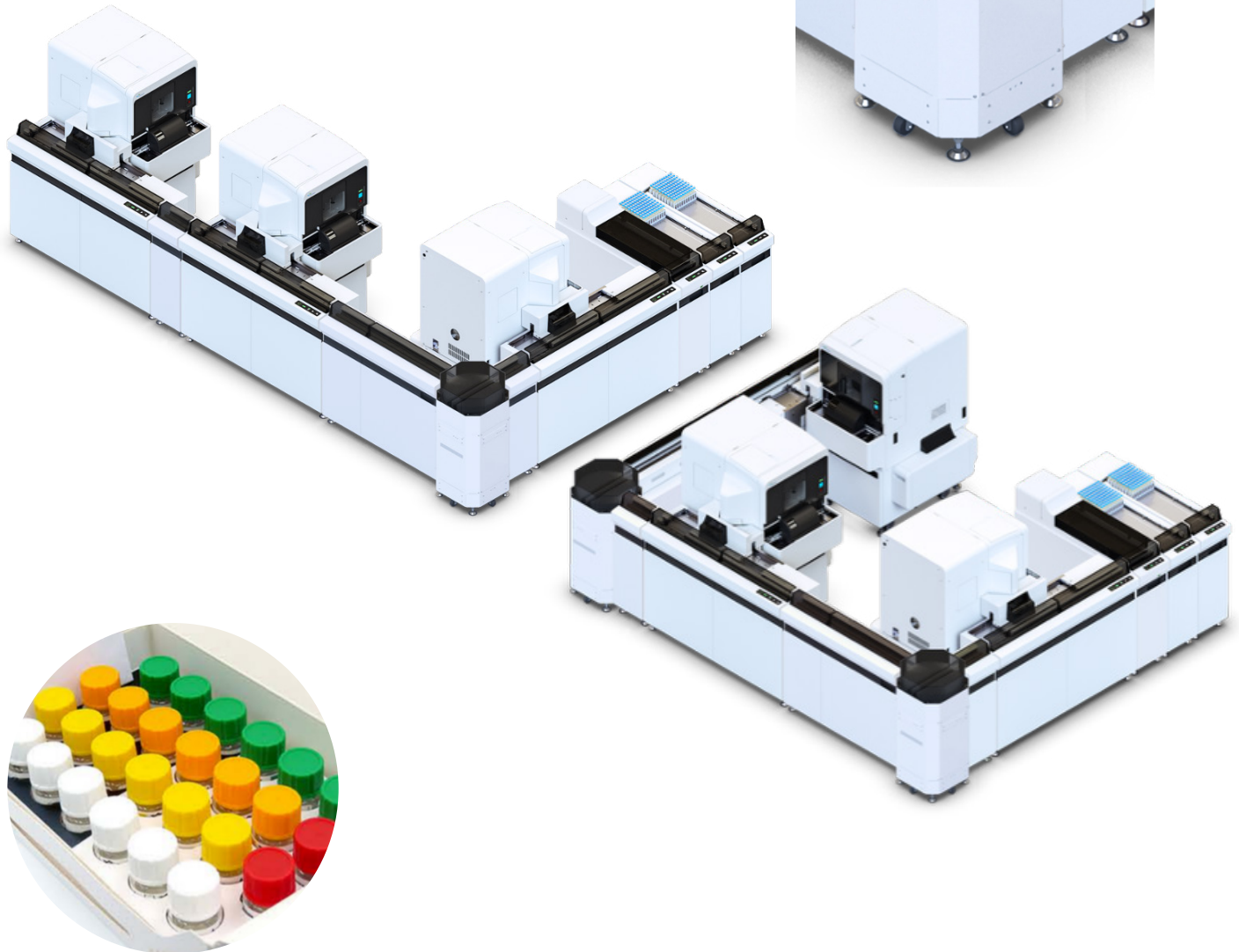

Industry tested and trusted reagents

Sysmex's hemostasis reagents and assay kits help labs like yours find the right solution for your patient testing. Providing proven reagents for routine testing, anticoagulant monitoring, bleeding risk assessment and thrombophilia screening, hemostasis solutions by Sysmex support screening, diagnosis and treatment of your most critical patients.

Contact your Sysmex Sales Representative to discover which reagents are right for your testing requirements.

Technical Specifications

The CN-9000 by Sysmex offers laboratories efficient workflows with the benefit of high-throughput and hands-free, walkaway operation. Connecting the CN-9000 to total lab automation can increase productivity in your hemostasis department, allowing staff more flexibility and time to perform other important tasks.

Analyzer Options	Maximum of three CN-6000 or CN-3000 analyzers TS-10 Automated Tube Sorter Turn Unit module for 90° left or right turn
Throughput	Process up to 1,000 samples/hour with three analysis modules
Data Storage	10,000 results/analyzer
Sampling	Continuous feed Cap piercing Dedicated STAT position Automated rerun and reflex testing

Sysmex Corporation

1-5-1-Wakinohama-Kaigandori Chuo-ku, Kobe 651-0073, Japan · Phone +81 78 265-0521 · www.sysmex.co.jp

Sysmex America, Inc.

577 Aptakisic Road, Lincolnshire, IL 60069, U.S.A. · Phone +1 888 879-7639 · www.sysmex.com/us

Sysmex Canada, Inc.

2680 Skymark Ave, Suite 1001, Mississauga, ON L4W 5L6 Canada · Phone +1 905 366-7900 · www.sysmex.ca

Sysmex Latin America and the Caribbean

Rua Joaquim Nabuco 615 - Bairro Cidade Jardim, São José dos Pinhais Paraná - Brasil - CEP 83040-210 · Phone +55 41 2104-1314 · www.sysmex.com.br